

Florence®

Florence is a chemically activated in office whitening gel containing 35% hydrogen peroxide. The new pH activated technology enables a fresh fast acting mix every time. This unique technology can be used to whiten teeth in a relatively short time. The system employs a novel Blockdam partially adhesive resin isolation technology to isolate the gingiva from the harmful effects of the peroxide ions. The whitening system does not require any light activation, making it the apt technology to be incorporated into every day dental practice.

Indications

- ◆ In office whitening for vital teeth.

Benefits

- ◆ No light activation.
- ◆ Relatively short application time.
- ◆ Novel isolation technology.
- ◆ Stabilized Hydrogen peroxide for better results and long shelf life.

- ◆ Easy application technique.

Chemistry

Hydrogen peroxide, alkali base.

Technical Data

1) PH: 9 - 10 %


Florence Bleaching Kit



Initial shade determination



**Application of Liquid dam.
Light cure for 20 seconds**

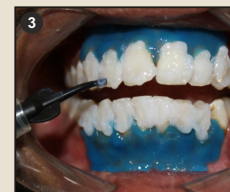

Application of Bleaching agent



**Shade change after 20 minutes
of bleaching**

Presentation

Ref. : 90001

Single Patient Kit :

- ◆ 2 x 0.3g Powder Pots
- ◆ 2 x 3ml Liquid Bottles
- ◆ 1 x 1.2ml Gingiva Shield VLC
(Gingival barrier Syringes)
- ◆ 2 x Micro Brushes
- ◆ 2 x Dispensing Tips
- ◆ 1 x Instruction Manual

Ref. : 90002

Two Patient Kit :

- ◆ 4 x 0.3g Powder Pots
- ◆ 4 x 3ml Liquid Bottles
- ◆ 2 x 1.2ml Gingiva Shield VLC (Gingival Barrier Syringes)
- ◆ 1 x Mixing Pot
- ◆ 4 x Micro Brushes
- ◆ 4 x Dispensing Tips
- ◆ 1 x Instruction Manual

PREVEST DENPRO LIMITED

Export Promotion Industrial Park (EPIP), Bari Brahmana, Jammu-181133, India.

E-mail : customercare@prevestdenpro.com **Website :** www.prevestdenpro.com

Toll Free No. : 1800-180-7172 Tel.: 91 1923 222774/222779