

Fusion-i-Seal®

Light Curing Glass Ionomer Composite Cement

Fusion-i-Seal® is ready to use, one component, light curing, fluoride releasing, radiopaque glass ionomer cement for base/lining for amalgam, ceramic and composite restorations. It is also indicated for extended fissure sealing.



Indications

- Lining under all kinds of filling materials
- Extended fissure sealing
- Restoration of smaller lesions

Benefits

- Long term fluoride release minimizes secondary caries
- Time saving application - Ready to use paste
- Prevents formation of marginal gaps thus protect from bacteria penetrations
- Protects the pulp
- Suitable for temporary and minor fillings
- Higher compressive and transverse strength surpasses ISO Standards

Chemistry

Fusion-I-Seal is fluoride ions containing ionomer glass in photo - polymerizable resin-difunctional monomers matrix.


Technical Data

- 1) Curing Time : 30 - 40 sec.
- 2) Flexural Strength : 80 MPa.

Application Technique



Cavity Preparation

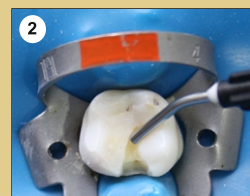

Placement of iséal



Light curing for 30 sec.



Completed liner placement

Presentation

- Ref. : 30002-1 (Trial Pack : 1 x 2g Syringe, 5 x Applicator Tips.)
 Ref. : 30002 (Intro Pack : 2 x 2g Syringes, 10 x Applicator Tips.)
 Ref. : 30003 (Economy Pack : 4 x 2g Syringes, 20 x Applicator Tips.)
 Ref. : 30003-1 (Intro Pack : 1 x 3g Syringe, 10 x Applicator Tips.)
 Ref. : 30003-2 (Kit Pack : 2 x 3g Syringes, 20 x Applicator Tips.)

PREVEST DENPRO LIMITED

Export Promotion Industrial Park (EPIP),
Bari Brahmana, Jammu-181133, India.

Tel.: 91 191 2430442/ 2432548
 Tel.: 91 1923 222774/222779
 Toll Free No. : 1800-180-7172

E-mail : info@prevestdenpro.com
 Website : www.prevestdenpro.com
 E-mail : customer@prevestdenpro.com

Fusion-i-Seal®

Light Curing Glass Ionomer Composite Cement



Fusion-i-Seal® is ready to use, one component, light curing, fluoride releasing, radiopaque, glass ionomer cement for base/lining for amalgam, ceramic and composite restorations. It is also indicated for extended fissure sealing.

PREVEST DENPRO LIMITED

Export Promotion Industrial Park (EPIP),
Bari Brahmana, Jammu-181133, India.

Tel.: 91 191 2430442/ 2432548
Tel.: 91 1923 222774/222779
Toll Free No. : 1800-180-7172

E-mail : info@prevestdenpro.com
Website : www.prevestdenpro.com
E-mail : customer@prevestdenpro.com