

# Orafil-G®

Temporary Filling Material, Paste

**Orafil-G** is a self curing, ready to use, radiopaque cement for provisional filling of dental cavities. Orafil G is a zinc sulphate - zinc oxide based cement and is designed for short term temporary applications.

This product is also available with Sodium Fluoride as Orafil Plus



## Indications

- Temporary restoration of cavities
- Temporary restoration between multiple endodontic visits
- Temporary restoration for inlay prepared cavity

## Benefits

- Non eugenol formulation offers non irritating properties
- Resin addition gives high adhesive properties
- Quick and void free curing in moist environment
- Pulp protecting and bactericidal
- Slight expansion of filling material ensures well sealed margin
- Excellent adaptability
- Simple to apply with filling instrument


## Chemistry

Zinc oxide, Zinc sulphate in EVA Resin/Glycol base with spearmint aroma

## Technical Data

- 1) Flow : 16 mm
- 2) Setting Time : 6min. 30sec.
- 3) Compressive Strength : 7-8 MPa


## Application Technique

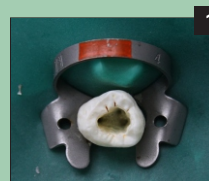


Endodontic Access



Placement of Material

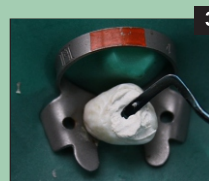

Material Adaptation

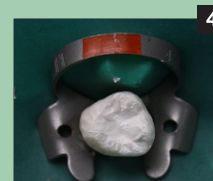

Completed Temporary Restoration

## Presentation

### Orafil G

Ref. : 11001  
1 x Jar 40g, 1 x Instruction Manual

### Orafil Plus

Ref. : 11002  
1 x Jar 40g, 1 x Instruction Manual

## PREVEST DENPRO LIMITED

Unit II, Export Promotion Industrial Park (EPIP),  
Bari Brahmana, Jammu-181133, India.

Tel.: + 91 191 350 6858, 60, 61  
Tel.: + 91 8899074151

E-mail : [info@prevestdenpro.com](mailto:info@prevestdenpro.com)  
Website : [www.prevestdenpro.com](http://www.prevestdenpro.com)

# Orafil-G®

Temporary Filling Material, Paste



**Orafil-G** is a self curing, ready to use, radiopaque cement for provisional filling of dental cavities. Orafil G is a zinc sulphate - zinc oxide based cement and is designed for short term temporary applications.

## PREVEST DENPRO LIMITED

Unit II, Export Promotion Industrial Park (EPIP),  
Bari Brahmana, Jammu-181133, India.

Tel.: + 91 191 350 6858, 60, 61  
Tel.: + 91 8899074151

E-mail : [info@prevestdenpro.com](mailto:info@prevestdenpro.com)  
Website : [www.prevestdenpro.com](http://www.prevestdenpro.com)